Flowlight LED Mask Pro

FLOWLIFE

Brochure

Flowlight LED Mask Pro

Red Light Therapy Mask

Table of Content

4 —	Manifest
6 —	Why Flowlight LED Mask Pro
7 —	Accessories
8 —	Flowlight LED Mask Pro - A Comprehensive Guide
10 —	Starting your Journey
16 —	Sleep Quality & Circadian Support
20 —	Stress Relief & Nervous System Reset
22 —	Energy & Cognitive Clarity
24 —	Skin Health & Facial Recovery
26 —	Contact Information
28 —	Manual
	• English • Svenska • Danska • Norsk • Suomi • Deutsch • Polski • Français • Italiano • Español • Nederlands

Manifest

In a world where limits are meant to be shattered, we fuel the fire of relentless pioneers. Our brand champions those who break the mold and embrace the extraordinary, symbolizing a life lived without compromise.

Committed to relentless innovation, we design for those who demand excellence and performance, understanding that the journey into the unknown often leads to unforgettable achievements.

We want to power dreams, aspirations, and triumphs of the audacious, sparking the passions of the unstoppable.

For those who redefine boundaries and explore new horizons, we are Flowlife. Made Possible.

Why Flowlight LED Mask Pro?

- **Targeted Light Therapy for Systemic Recovery** - Flowlight LED Mask Pro is a wearable red light therapy device delivering 630 nm and 850 nm wavelengths—clinically proven to stimulate mitochondrial activity, reduce inflammation, and accelerate recovery. By concentrating light on the most vascular and neurologically rich surface in the body—the face—it promotes cellular repair throughout the nervous system, improves circulatory balance, and helps regulate hormonal and autonomic functions—a local application with full-body impact.
- **Full-Face LED Architecture with Clinical Efficiency** - Each of the 72 LEDs emits both wavelengths through a wide 120° beam angle, ensuring complete and even coverage across the face. This 3-in-1 diode design maximizes light saturation in every zone, from forehead to jawline, improving penetration and energy absorption without gaps. With a controlled irradiance of ~30 mW/cm², it delivers a therapeutically effective dose—strong enough to trigger a biological response, gentle enough for consistent daily use.
- **Three Modes, Three Timers, Three Intensities** - The wireless controller offers precision without complexity. Choose between combined red + NIR, red only, or NIR only, depending on your recovery goal—whether it's parasympathetic activation, skin-level repair, or deeper tissue regulation. Each mode offers three-time settings (5, 10 or 15 minutes) and three intensity levels, making the Flowlight LED Mask Pro adaptable to your energy levels, schedule, and desired depth of stimulation.
- **Portable Recovery, Built for Active Lifestyles** - Flowlight LED Mask Pro isn't just for the skin—it's a systemic recovery tool designed for daily use and full-body impact. While it doesn't replace full-panel therapy, it complements it by targeting the most neurologically sensitive zone in the body. Use it as a standalone tool or integrate it into your broader recovery ecosystem—at home, on the move, or between sessions. It's clinical-grade recovery, simplified, and made wearable.


Accessories

Flowlight LED Mask Pro comes with everything you need for consistent, targeted light therapy, whether at home or on the go. Each item is selected to support your treatment routine and keep the device protected, charged, and ready to use.

Flowlight LED Mask Pro

The core device. Designed with wide-angle 3-in-1 LEDs delivering both 630 nm and 850 nm wavelengths. It offers complete control via the included wireless controller or simplified access through its built-in power button for default mode—lightweight and contour-fit for facial ergonomics.

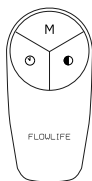Wireless Controller

Optional for users who want full treatment customization. Adjusts mode, time, and intensity with responsive tactile buttons and visual indicators for confirmation of settings. Enables three light programs, three-time durations (5/10/15 minutes), and three intensity levels.


USB-C Charging Cable

Use to charge your Flowmask.

Protective Carry Pouch

Slim, durable, and structured for storage or travel. Keeps the mask and accessories protected without bulk.

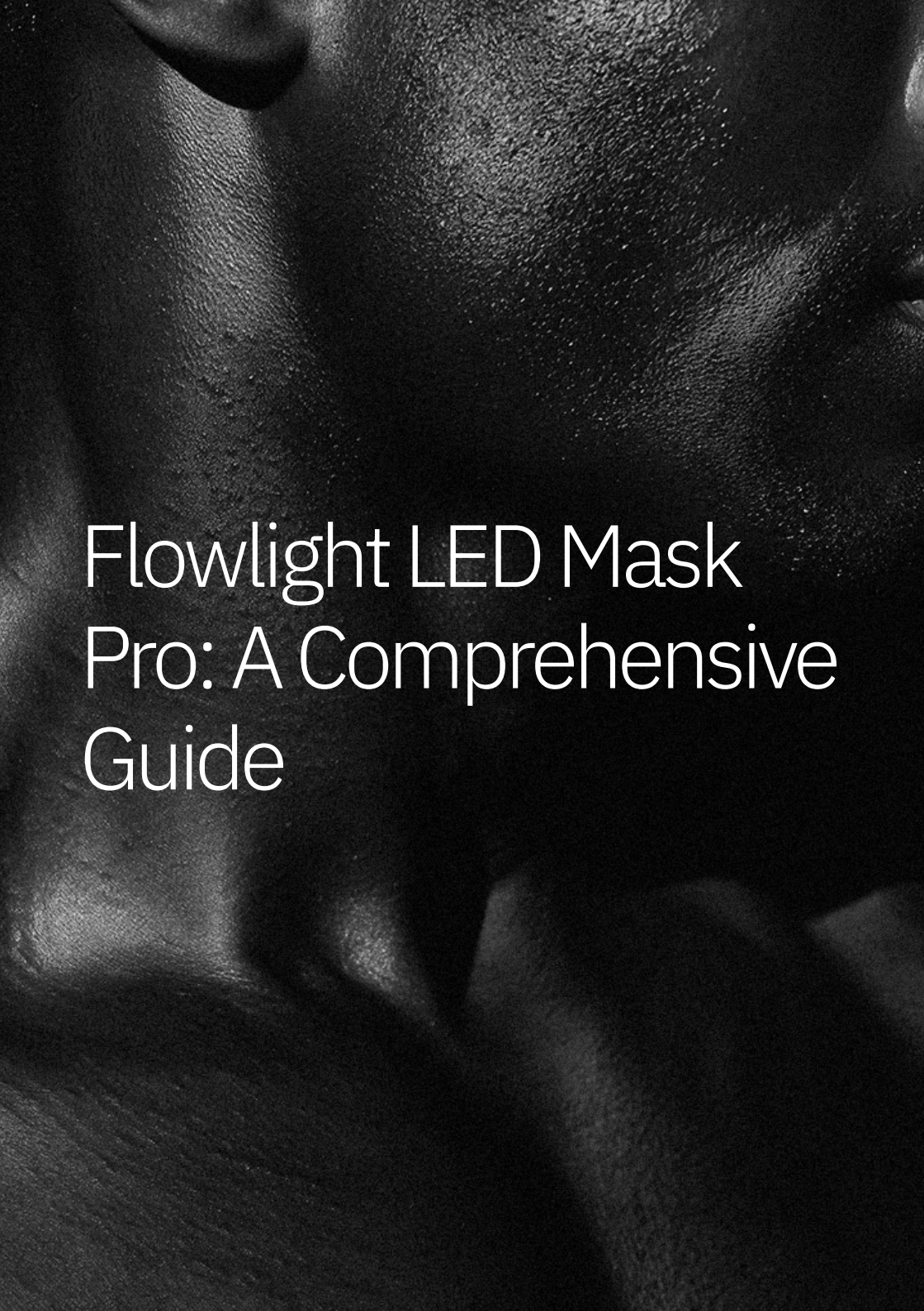A high-contrast, black and white close-up photograph of a person's face. The lighting is dramatic, coming from the side, which highlights the texture of the skin and the contours of the face. The person's eyes are closed, and their mouth is slightly open. The overall mood is artistic and moody.

Flowlight LED Mask Pro: A Comprehensive Guide

To make the most of your sessions, it's essential to understand how Flowlight LED Mask Pro works and how to use it effectively, both with and without the wireless controller. This guide walks you through the entire experience, from activation to storage.

Preparing for Your Session

Start with clean, dry skin

Ensure your face is clean and free from make-up, oil, and lotion. Clean skin maximizes light absorption and minimizes interference with the mask's optical output.

No need for eye protection

Flowlight LED Mask Pro is designed with a safe, clinically aligned wavelength distribution. Most users prefer to close their eyes during treatment for comfort, but protective glasses are not required.

Starting Your Session

Activate Default Mode (manual use)

Hold the small button located on the top-right inner side (above your right eye when worn) for 3 seconds. This activates the device and the default mode.

- 15-minute treatment.
- Combined 630 nm red and 850 nm near-infrared (NIR wavelengths).
- Intensity level 3 (highest).

Note: Without the controller, only the Default Mode is available. You cannot adjust mode, intensity, or time.

Understanding the Modes

To customize your session, connect the wireless controller. This unlocks complete control of the following:

- Mode: Red + NIR / Red only / NIR only.
- Timer: 5, 10, or 15 minutes.
- Intensity: Level 1 (Low), 2 (Medium), or 3 (High).

1. Red + NIR (Default Mode)

Wavelengths: 630 nm, 850 nm

Delivers broad-spectrum stimulation from superficial dermal layers to deeper neurological tissues. Supports systemic recovery, mitochondrial repair, and the regulation of inflammation.

2. Red Only (630 nm)

Optimized for skin-level applications. Stimulates collagen production, accelerates wound healing, and supports mood regulation.

3. Near-Infrared Only (850 nm)

Penetrates deeper into the facial nerves, vasculature, and connective tissues. Ideal for nervous system regulation, hormonal balance, and fatigue relief.

Intensity Levels

Each mode can be adjusted to 3 intensity levels via the controller. Higher levels increase irradiance but do not change wavelength.

Timer Settings

Available durations are 5, 10, or 15 minutes. The mask powers off automatically at the end of the session.

Visual & Sensory Mode Indicators

To help confirm which light mode is active, the Flowlight LED Mask Pro uses a built-in indicator light positioned on the outside of the mask's right-side module:

- White light = Mode 1 (Red + NIR)
- Solid red light = Mode 2 (Red only)
- Flashing red light = Mode 3 (NIR only)

In addition to the external indicator, the mask offers sensory and visual cues during treatment:

- Red-only mode emits a visible red light.
- Red + NIR mode emits the same red light, often perceived as slightly deeper due to the combined wavelengths.
- NIR-only mode may appear inactive to the eye, as near-infrared light is invisible. However, the mask will still emit gentle warmth and remain fully active. Use the indicator light to confirm mode selection.

Using the Wireless Controller

The Flowlight LED Mask Pro controller allows for complete customization of your treatment. It is wireless, intuitive, and syncs directly with the mask using both audible and visual feedback.

Initial Activation

After inserting the battery into the controller, hold the mode button "M" to activate the remote control.

Button Layout & Functions

Mode Button (top middle)

- Press once to cycle between light modes:
 - Mode 1: Red + NIR
 - Mode 2: Red only
 - Mode 3: NIR only
- The selected mode is confirmed by the indicator light on the outer surface of the mask's right-side module.

Note: This button does not power the mask on or off—it only changes the light mode.

Time Button (bottom left)

- Press to select treatment duration: 5, 10, or 15 minutes.
- The LED display shows 1, 2, or 3 lights, corresponding to the selected duration.

Intensity Button (bottom right)

- Press to set intensity level: Low (1), Medium (2), or High (3).
- The LED display shows 1, 2, or 3 lights to reflect the selected brightness.

Audible Feedback

Every adjustment—mode, timer, or intensity—is confirmed by a soft beep from the mask, ensuring real-time responsiveness even when the visual indicator is not in view.

During Your Session

Relax and disconnect

Flowlight LED Mask Pro was designed for portability, making it perfect for modern recovery routines on the go. But while you can technically use the mask on the go, we strongly recommend making each session intentional. Sit or lie down. Let go of stimulation.

Monitor comfort

If using the controller, you can adjust intensity or switch modes mid-session. The mask will continue to operate at the new setting until the timer runs out or you manually end the session.

After Your Session

Automatic shut off

The mask turns off automatically when the session ends. You can also stop it early by holding the power button on the mask.

Cleaning and storage

Gently wipe the device with a dry microfiber cloth. Store it in the provided carry pouch to keep it clean and safe.

Charging

Use the included USB-C charger to recharge the mask. Use two AAA-Battery for the remote control. One full charge typically supports multiple sessions.

Sleep Quality & Circadian Support

Optimizing both ends of the sleep cycle—
morning activation and evening regulation.

Morning: Set Your Internal Clock

Why it matters

Your circadian rhythm is primarily regulated by light. Early exposure to red light tells your brain it's time to be awake, setting the internal clock for melatonin production 12–14 hours later. This improves your ability to fall asleep naturally in the evening.

How it works

- Red light (630 nm) suppresses melatonin and elevates natural cortisol to promote wakefulness.
- Enhances morning mood, energy, and mental clarity.
- Strengthens circadian timing for more stable nighttime regulation.

Recommended Settings

Mode

- Red Only or Combined Mode (Red + Near-Infrared).

Intensity

- Level 2–3.

Duration

- 5–10 minutes

Time of Day

- Within 1 hour of waking.

Frequency

- Daily, especially in darker seasons or after travel.

Evening: Wind Down for Restorative Sleep

Why it matters

As evening approaches, your nervous system should gradually shift from alert to parasympathetic mode. Near-infrared light supports this transition by alleviating cranial tension and signaling the release of melatonin.

How it works

- Near-infrared light (850 nm) supports parasympathetic activation and reduces systemic stress.
- Calms facial and cranial muscles, often linked to poor sleep and headaches.
- Encourages melatonin release by avoiding light wavelengths that suppress it.

Recommended Settings

Mode

- Near-Infrared Only (for deep relaxation).

Intensity

- Start at Level 1–2; increase to Level 3 if well-tolerated.

Duration


- 10–15 minutes


Time of Day

- 30–60 minutes before bed.

Frequency

- Nightly or after high-stress evenings.

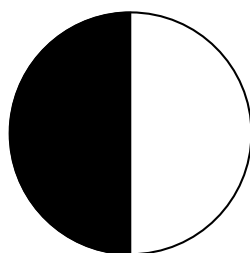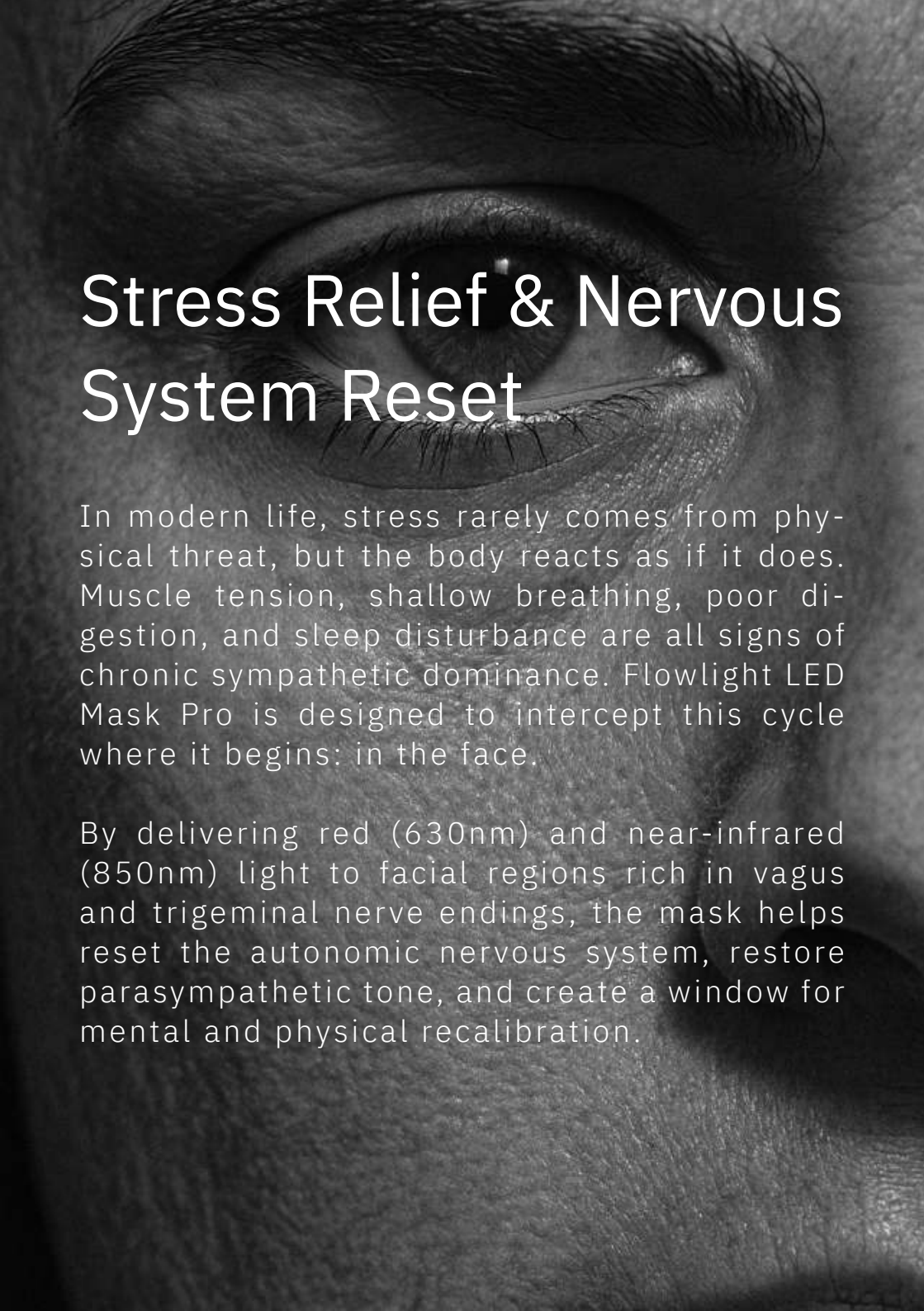Stress Relief & Nervous System Reset

In modern life, stress rarely comes from physical threat, but the body reacts as if it does. Muscle tension, shallow breathing, poor digestion, and sleep disturbance are all signs of chronic sympathetic dominance. Flowlight LED Mask Pro is designed to intercept this cycle where it begins: in the face.

By delivering red (630nm) and near-infrared (850nm) light to facial regions rich in vagus and trigeminal nerve endings, the mask helps reset the autonomic nervous system, restore parasympathetic tone, and create a window for mental and physical recalibration.

How it works

- Light absorption through cranial nerve branches (especially supraorbital and infraorbital) activates parasympathetic pathways.
- This slows the heart rate, lowers blood pressure, and reduces the levels of circulating stress hormones.
- Regular exposure may improve vagal tone, heart rate variability, and emotional regulation.
- Many users report a sense of calm afterward, often accompanied by improved sleep or digestive ease.

Recommended Settings

Mode

- Red + NIR (for dual action on surface and subdermal nerves)

Intensity

- Level 1–2 (gentle input is more effective for nervous system regulation).

Duration

- 10–15 minutes

Time of Day

- Afternoon or early evening—avoid use during multitasking or high-activity.

Frequency

- 3–7 times per week, depending on baseline stress levels.

Energy & Cognitive Clarity

Brain fog, low motivation, and afternoon crashes often stem from disrupted cellular metabolism and poor alignment with the circadian rhythm. The face—home to dense vascular networks and cranial nerve branches—provides a unique entry point for influencing systemic energy regulation. Flowlight LED Mask Pro channels clinically studied red and near-infrared light to stimulate mitochondrial activity, increase oxygen delivery, and elevate mental sharpness.

How it works

- Red and NIR light penetrates superficial vascular beds and cranial nerve endings, stimulating ATP production in facial tissues and adjacent areas of the frontal cortex.
- Enhanced mitochondrial output improves cellular energy availability, not just locally but systemically.
- Cortisol levels are modulated through hypothalamic pathways, promoting alertness without overstimulation.
- Users report clearer thinking, improved focus, and stabilized energy throughout the day.

Recommended Settings

Mode

- Red only or Red + NIR (for surface activation and deeper circulation support).

Intensity

- Level 2–3 (depending on sensitivity and desired stimulation)

Duration

- 5–15 minutes.

Time of Day

- Morning to early afternoon—avoid using late in the day to prevent interference with sleep.

Frequency

- 3–6 times per week, ideally paired with a morning ritual or high-focus block.

Skin Health & Facial Recovery

Flowlight LED Mask Pro delivers clinically validated wavelengths that support collagen production, reduce inflammation, and restore skin vitality. Unlike topical treatments, red and near-infrared light reaches deep into the dermis and subcutaneous layers, stimulating cellular renewal from within..

How it works

- Red light (630 nm) activates fibroblasts to boost collagen and elastin production.
- Near-infrared light (850 nm) increases microcirculation and cellular metabolism in deeper tissues.
- Anti-inflammatory effects calm reactive skin and accelerate recovery from irritation, acne, or fatigue.
- Repeated exposure helps maintain skin firmness, balance hydration, and achieve an even tone over time.

Recommended Settings

Mode

- Red only (surface repair) or Red + NIR (for combined dermal and subdermal support).

Intensity

- Level 2 for daily use
- Level 3 for deeper treatment or post-stressor repair.

Duration

- 5-10 minutes, depending on skin sensitivity.

Time of Day

- Morning or post-cleansing routines; avoid immediately before intense sun exposure.

Frequency


- 5–7 times per week for active improvement
- 3–4 times per week for maintenance

hello@flowlife.com
www.flowlife.com
Luntmakargatan 90
113 51 Stockholm

Manual

Flowlight LED Mask Pro

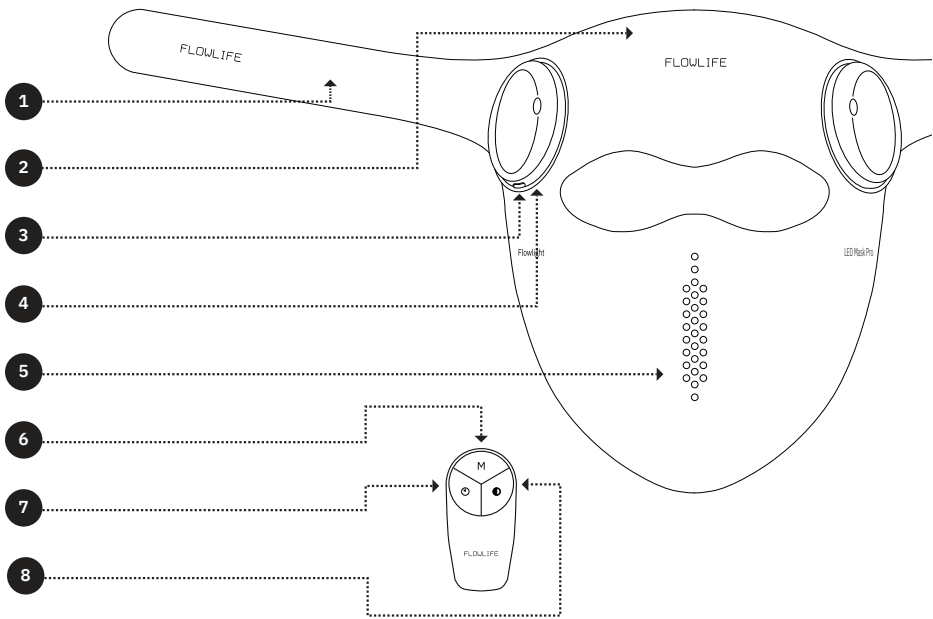1. Velcro Strap
2. Mask
3. On/Off
4. USB Charging Port
5. Breathing Vents
6. Mode Selection
7. Timer
8. Light Intensity

ENGLISH SAFETY INSTRUCTIONS

ELECTRICAL HAZARDS WARNINGS

- Never clean the mask while plugged in
- Do not use near water (bath, sink, etc.)
- Unplug immediately if it comes into contact with water

USAGE WARNINGS

- Never cover the mask during use
- Do not use on broken or irritated skin
- Only use original accessories
- Keep out of reach of children
- Never insert objects or attempt repair
- Store away from hot surfaces or damp areas
- Do not use while sleeping, driving, or operating machinery

STOP USE IF

- The product or cable is damaged
- It has been dropped or exposed to water

USAGE

- Ensure your face is clean and dry before use
- Place the mask securely on your face for full contact
- Fasten the Velcro strap behind your head to keep it in place during use
- Press and hold the small button on the inside top-right module (above your right eye when worn) for 3 seconds to activate the mask

DEFAULT MODE (MANUAL START)

- Duration: 15 minutes
- Wavelengths: 630 nm red + 850 nm near-infrared
- Intensity: Level 3 (highest)

WIRELESS CONTROLLER (OPTIONAL FOR FULL CUSTOMIZATION)

- Insert 2 × AAA batteries into the wireless controller before first use
- After installing the batteries, all 6 indicator lights on the controller will flash 5 times—this confirms successful power-up and connectivity

MODE BUTTON (TOP MIDDLE)

- Press once to cycle between the three light modes:
- → Mode 1: Red + NIR
- → Mode 2: Red only
- → Mode 3: NIR only
- Mode indicator (on right-side module of mask):
- – White light = Mode 1 (Red + NIR)
- – Solid red = Mode 2 (Red only)
- – Flashing red = Mode 3 (NIR only)
- Audible beep confirms each selection
- Note: This button does not power the mask on or off

TIME BUTTON (BOTTOM LEFT)

- Press to set session duration: 5, 10, or 15 minutes
- LED indicators:
- – 1 light = 5 minutes
- – 2 lights = 10 minutes
- – 3 lights = 15 minutes

INTENSITY BUTTON (BOTTOM RIGHT)

- Press to set intensity level: Low, Medium, or High
- LED indicators:
- – 1 light = Low
- – 2 lights = Medium
- – 3 lights = High

DURING USE

- Sit or lie down comfortably
- Avoid multitasking
- Keep eyes closed if preferred

AFTER USE

- The mask will shut off automatically when the session ends
- To end manually, hold the button on the mask

CLEANING

- Wipe both inner and outer surfaces with a dry microfiber cloth after each use

CHARGING

- Use the included USB-C cable to charge the mask
- Charging indicator (on mask):
 - Red flashing light = Charging
 - White light = Fully charged

RISK GROUPS – ALWAYS CONSULT A DOCTOR IF:

- You have had recent burns, fever, infections, or photosensitivity
- You are pregnant (avoid the abdominal area)
- You are diagnosed with SLE, cancer, hyperthyroidism, or severe bleeding disorders
- You use photosensitizing medications
- You have eye disease or visual sensitivity
- You have a history of melasma or hyperpigmentation

DISCLAIMER

- This product is not a medical device
- It is designed for non-medical home and professional wellness applications
- Effects vary between individuals
- Consult a healthcare provider for specific concerns

ELECTRONIC WASTE & DISPOSAL

- This product contains a lithium-ion battery and electronic components
- Do not dispose of in household waste
- Recycle electronic waste in accordance with local regulations

WARRANTY

- A 2-year limited warranty, effective from the date of purchase, covers the Flowlight LED Mask Pro
- Covers defects in materials and workmanship under regular use
- Does not cover damage resulting from misuse, unauthorized repairs, accidents, or exposure to water
- Proof of purchase is required for warranty claims

TECHNICAL SPECIFICATIONS

- Wavelengths: 630nm (LED RED), 850nm (LED NIR)
- Optical Power per Diode:
 - 630nm: 6mW
 - 850nm: 15mW
- Beam Angle: 120°
- Irradiance on Skin: ~30 mW/cm²
- Battery: 1440mAh lithium-ion (40 min operation per charge)
- Number of Diodes: 216 pcs (144 x 630nm + 72 x 850nm)

OPERATION MODES

- Mode 1: Red + NIR combined
- Mode 2: Red only
- Mode 3: NIR only

TIMER OPTIONS

- 5, 10, or 15 minutes

INTENSITY LEVELS

- Low, Medium, High

CHARGING INDICATOR LIGHT

- The indicator will flash red when charging. And full white light when fully charged.

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
USB-C OPLAADKABEL

- Nominale ingang: AC100–240V
- Nominale uitgang: DC 5V–1A

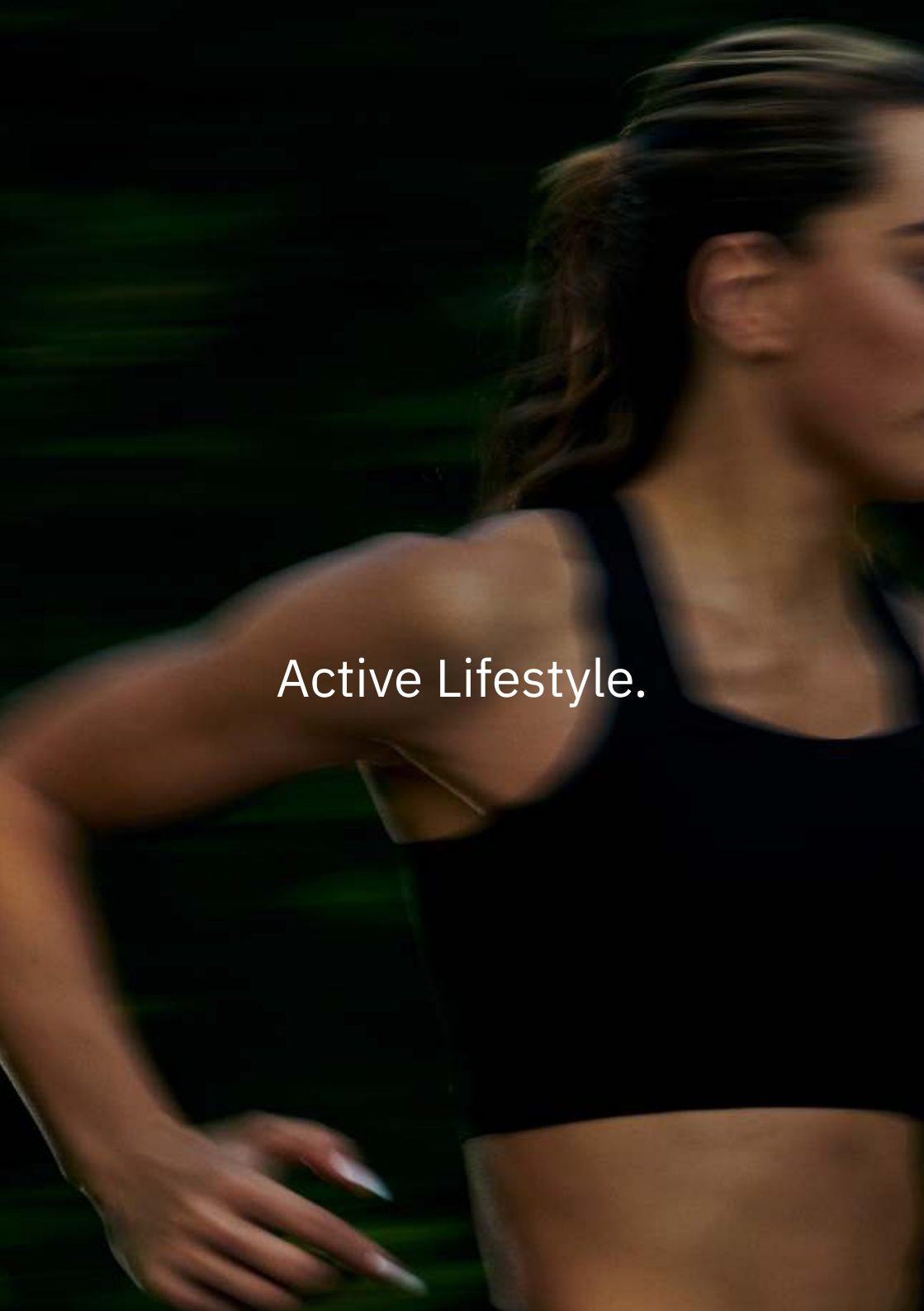A close-up, side-profile shot of a woman running at night. She is wearing a black athletic top. Her hair is dark and pulled back, with some strands blowing in the wind. The background is dark and blurred, suggesting a nighttime outdoor setting. The text "Active Lifestyle." is overlaid in white, sans-serif font in the center of the image.

Active Lifestyle.

Made Possible.

A black and white photograph showing a close-up of a person's neck and shoulder. The person is wearing a dark, possibly black, strap over their shoulder. The skin is smooth and has a slight sheen. The word "FLOWLIFE" is overlaid in white, uppercase, sans-serif font across the center of the image.

FLOWLIFE