
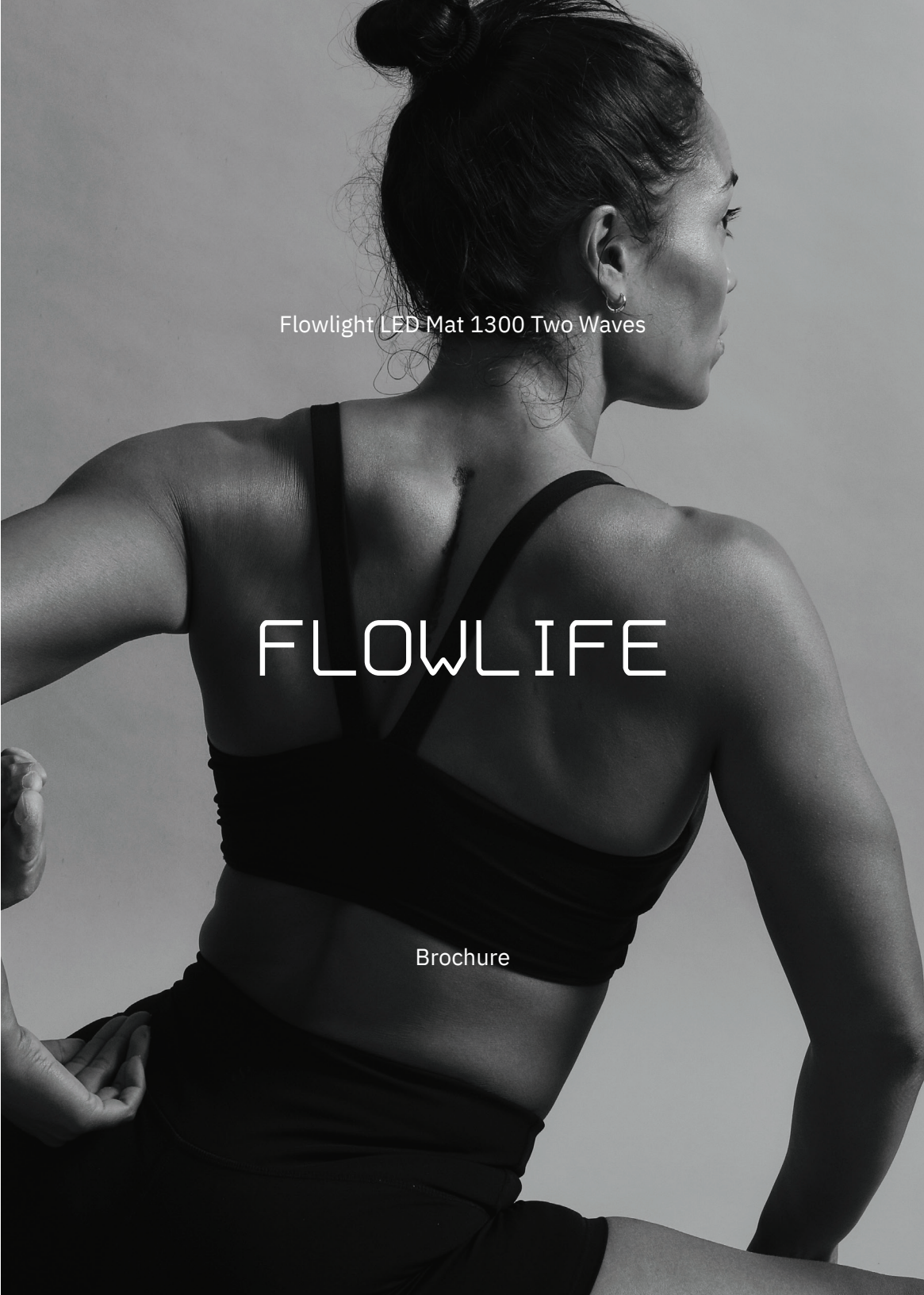
Flowlight LED Mat 1300 Two Waves

FLOWLIFE

Brochure

Flowlight LED Mat 1300 Two Waves

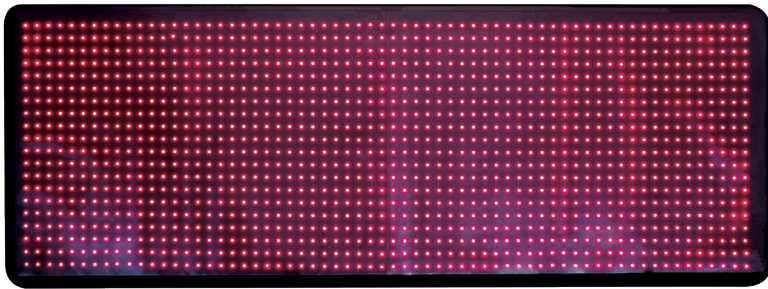Red Light Therapy Mat

Table of Content

4 —	Manifest
6 —	Why Flowlight LED Mat 1300 Two Waves
7 —	Accessories
8 —	Flowlight Mat – A Comprehensive Guide
10 —	Starting Your Journey Understanding the Light Modes
15 —	Using the Controller
18 —	Muscle Soreness and Recovery
22 —	Reducing Chronic Pain and Inflammation
28 —	Enhancing Skin Health and Cellular Repair
33 —	Contact Information
34 —	Manual • English • Svenska • Danska • Norsk • Deutsch • Français • Español • Italiano

Manifest

In a world where limits are meant to be shattered, we fuel the fire of relentless pioneers. Our brand champions those who break the mold and embrace the extraordinary, symbolizing a life lived without compromise.

Committed to relentless innovation, we design for those who demand excellence and performance, understanding that the journey into the unknown often leads to unforgettable achievements.

We want to power dreams, aspirations, and triumphs of the audacious, sparking the passions of the unstoppable.

For those who redefine boundaries and explore new horizons, we are Flowlife. Made Possible.

Why Flowlight LED Mat 1300 Two Waves Two Waves?

The Flowlight LED Mat 1300 Two Waves is a next-generation red light therapy mat designed to treat the whole body with high-density dual-wavelength LED light. With 3840 total LEDs across a flexible, ergonomic surface, it delivers therapeutic intensity while remaining soft and adaptable.

It's more than just a light mat. It's a foundation for everyday recovery.

High-Density Coverage

130 mW/cm² at skin level for deep-tissue penetration and high output per square centimeter.

Dual Wavelength Precision


Targets both superficial and deep tissue layers with clinically proven 660nm red and 850nm NIR light.

Full-Body Recovery Anywhere

Flexible, foldable, and lightweight. Designed for athletes, therapists, and anyone seeking serious recovery from home.

Accessories


Everything you need to get started is included.

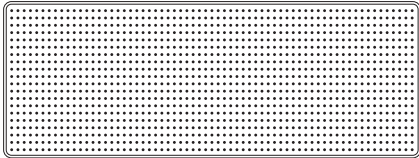
Flowlight LED Mat 1300
Two Waves LED Mat


Color-coded remote controller

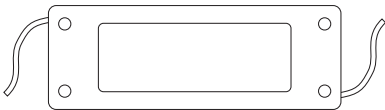UL-certified power driver (24V5A)

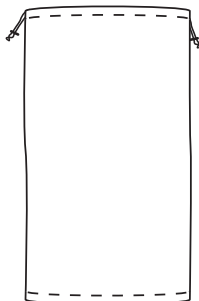
Protective glasses

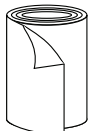Vetrol straps

Travel carry bag

Instruction manual

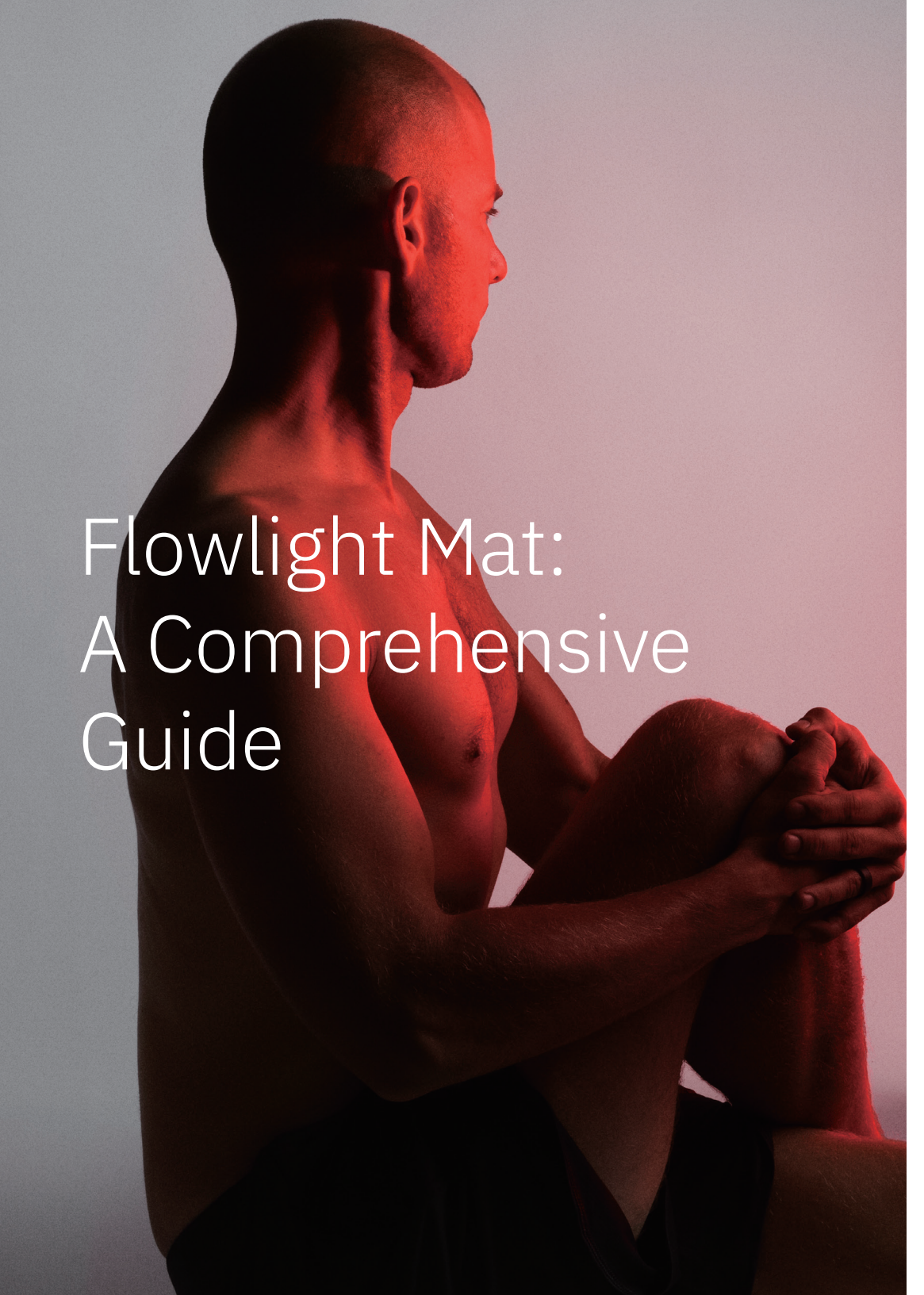A full-page background image featuring a man in a yoga pose, specifically a seated forward bend (Uttanasana). He is shown from the side, with his torso twisted towards his right leg. His hands are clasped together, holding his right foot. The lighting is a strong, warm red, creating a dramatic silhouette effect and highlighting the contours of his muscles. The background is a plain, light-colored wall.

Flowlight Mat: A Comprehensive Guide

This mat combines science and simplicity. With 3840 LED chips, a breathable textile shell, and full dimming and pulsing control, it's built to treat muscle soreness, joint stiffness, skin health, and sleep regulation through optimized photobiomodulation. Use it lying down on a massage table, or post-training as part of your recovery routine. You control every parameter: intensity, wavelength, pulsing frequency, and session duration.

Starting Your Journey

Prepare Your Body and Environment

Before your first session, take a few minutes to prepare both your body and your space. A calm, clean setup helps ensure the light reaches your skin effectively and supports both physical and mental recovery.

Recommendation

Use the mat on a firm, flat surface such as a massage table, yoga mat, or floor. Your skin should be clean, dry, and free from oils or lotions, as these can reduce light penetration. Remove thick clothing from the treatment area to allow for direct light exposure. If you're treating near the face, use the included protective glasses. Silence your phone, dim the room lights, and allow yourself to settle in.

Understanding Light Modes: RED, NIR & Combined

Flowlight LED Mat 1300 Two Waves combines cutting-edge red and near-infrared (NIR) light therapy for a comprehensive full-body wellness experience. By utilizing two therapeutic wavelengths, 660nm red light and 850nm near-infrared light, each mode delivers tailored benefits for skin regeneration, deep tissue recovery, and overall health.

RED/NIR Mode (Combined: 660nm + 850nm)

In **RED/NIR Mode**, both 660nm red light and 850nm near-infrared light are activated simultaneously. This combination is ideal for full-body recovery, as it addresses both surface and deeper tissues. Red light works at the surface level, while near-infrared light penetrates deep into muscles and joints, enhancing recovery.

660nm red light:

This wavelength improves skin health, promotes collagen production, and helps reduce surface inflammation. It stimulates the superficial tissues, increasing blood flow and oxygenation in the skin, improving tone, and rejuvenating the surface.

850nm near-infrared light:

It penetrates deeper into muscles and joints, stimulating cellular repair and promoting better circulation. It supports muscle recovery, reduces soreness, and accelerates tissue healing in deeper layers that red light cannot reach.

Recommended for: Full-body recovery

This mode is perfect for treating both superficial and deep tissue issues, making it suitable for a wide range of recovery needs.

Post-exercise recovery

It's ideal after a workout or other physically demanding activity, as it supports muscle repair and reduces soreness.

Systemic inflammation

If you're dealing with widespread inflammation, this combination of red and NIR light can provide relief.

Chronic pain management

It helps with chronic pain, including back pain, muscle stiffness, and arthritis.

Ideal time of day:

Morning or midday:

This mode delivers a comprehensive, energizing treatment ideal for starting the day or rejuvenating your body throughout the day. You can also use it in the evening at a lower intensity for deep relaxation and sleep support.

RED Mode (660nm only)

Red light at 660nm targets the upper layers of the body's tissue – primarily the epidermis and dermis, but its physiological effects go beyond skin-deep. It stimulates cellular function in surface tissue, enhances local blood flow, and modulates inflammation. This mode is particularly effective for targeting areas near the skin, such as joints, superficial muscles, and connective tissue.

660nm red light

is absorbed by the mitochondria in surface-level cells, increasing ATP (cellular energy) production and improving cell signaling. As a result, it promotes both visible and functional repair: from collagen regeneration to anti-inflammatory effects.

Recommended for:

Skin health and repair

Stimulates fibroblasts and supports collagen and elastin synthesis, reducing fine lines, uneven texture, and dryness.

Surface inflammation

Helps calm inflamed tissue from irritation, friction, or repetitive stress injuries.

Superficial muscle zones

Effective on areas like the neck, shoulders, forearms, and joints, where muscles are closer to the skin surface.

Microcirculation

Increases blood flow to capillaries, supporting faster nutrient delivery and waste removal in local tissues.

Tissue oxygenation

Enhances local oxygen exchange, helping sore areas recover more efficiently.

Ideal time of day:

Morning or early afternoon:

Red light can stimulate your skin and circulation, which may not be conducive to a relaxing evening. Use this mode earlier in the day when you need an energy boost and improved skin health.

Setting the Timer

Timer Range

The Flowlight LED Mat 1300 Two Waves allows you to set treatment duration anywhere from 1 to 30 minutes using the controller. The session time is displayed on the controller and counts down automatically.

Recommendation

For first-time users, start with 10–15 minutes per session. As your body adapts, you can gradually increase the session length to 20–30 minutes, depending on your treatment goal.

Short sessions (10–15 min):

Pre-workout activation or targeted use.

Medium sessions (15–25 min):

General recovery and stress relief.

Extended sessions (25–30 min):

Chronic pain, deep recovery, or skin-focused treatments.

NIR Mode (850nm only)

NIR Mode utilizes the 850nm wavelength of near-infrared light. Unlike visible red light, near-infrared light is invisible but penetrates deeper into tissues. This mode is designed for muscle recovery, joint pain relief, and deep-tissue healing. NIR light stimulates the mitochondria in cells, accelerating the repair and regeneration of tissue.

850nm near-infrared light:

This wavelength penetrates deeply, reaching muscle fibers, tendons, and joints. It helps to accelerate healing by improving cellular energy production (ATP), promoting better circulation, and reducing inflammation in deeper tissues. This makes it ideal for treating chronic pain and deep muscle soreness, and for promoting healing in larger muscle groups.

Recommended for:

Deep tissue recovery

This mode is ideal for muscle soreness, joint stiffness, and post-exercise recovery. It enhances oxygen and nutrient delivery to muscle tissue, accelerating repair.

Chronic pain and joint stiffness

It's especially beneficial for conditions such as arthritis, tendinitis, and other inflammatory conditions affecting the joints and muscles.

Muscle and nerve recovery

Deep penetration helps to repair muscle fibers and improve flexibility, while also soothing nerve pain.

Ideal time of day:

Evening:

Use NIR mode in the evening to help calm your nervous system, reduce cortisol levels, and promote better sleep. NIR light also enhances melatonin production after the session, supporting deeper sleep.

Post-workout:

Great for relieving soreness and improving recovery after intense physical activity.

Using the Controller


The Flowlight LED Mat 1300 Two Waves is designed to be intuitive and highly customizable. Its external controller gives you complete control over how each session is delivered. Whether you prefer minimal input or want to fine-tune every variable, the mat adapts to your needs.

Here's how to operate and optimize your treatment settings.

Powering On

1. Plug the mat into the power supply and connect the controller.
2. Connect the male connector of the remote cable to the cable attached to the Flowlight Mat. Make sure the slots on both connectors are aligned, then push them together firmly and secure the connection by tightening the locking ring.
3. Connect the female connector of the remote cable to the power adapter cable. Align the slots, press firmly until fully connected, and tighten the locking ring.
4. Plug the power adapter into a standard wall outlet.
5. When power cord is connected to a wall socket, the interface screen lights up. If you want to shut it down, unplug the power cord from the wall socket.
6. Choose your desired settings before starting the session.

Key Controls Overview

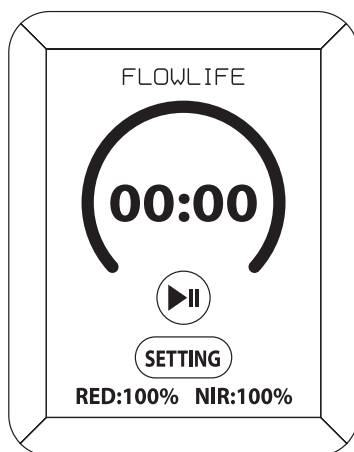Function Description

- Timer - Set the treatment duration: from 1 to 30 minutes
- Intensity - Adjust each wavelengths brightness from 1% to 100% , depending on comfort and goal
- Pulse - Set pulse frequency (1–1000Hz) for rhythmic on/off light modulation

Setting the Timer

Timer Range

The timer can be set for 1 to 30 minutes, depending on your treatment needs. The countdown is displayed on the controller, and the mat automatically turns off when the session ends.

Adjusting Light Intensity

Light Intensity

The dimmer adjusts the light brightness from 1% to 100%. Higher intensity delivers deeper penetration and greater therapeutic effect.

Recommendation

Use 100% intensity for maximum results in most cases. If you're new to red light therapy, treating sensitive areas, or using the RED or RED/NIR mode before sleep, start with 50–80%. Increase gradually as your body adjusts.

Pulse Setting

Pulse Range

The mat allows you to set pulsing frequencies between 1Hz and 1000Hz. Pulsing modulates light to turn it on and off rhythmically, which may enhance cellular responses and tissue stimulation.

Recommendation

Use low pulse settings (1–10Hz) to promote relaxation and nervous system regulation. For muscle recovery, moderate to high frequencies (40–300Hz) may increase circulation and support more profound healing. Higher settings (500–1000Hz) are recommended for advanced use and can enhance tissue penetration.

During Your Session

Relax and Recover

The Flowlight LED Mat 1300 Two Waves is designed for passive, full-body recovery. Once activated, lie back and allow the light to reach your skin without barriers. Maintain direct contact with the mat for optimal light absorption.

Recommendation

Wear lightweight or no clothing on the treatment area. Keep your breathing steady and avoid distractions. Dim surrounding lights and allow your body to relax fully. You may feel gentle warmth, but the mat should never feel uncomfortably hot. Adjust intensity if needed.

Ending Your Session

Power Down Safely

At the end of your session, the mat turns off automatically if the timer has been used. If not, manually power off the device using the controller before disconnecting the power supply.

Recommendation

Let the mat cool for a few minutes before handling. Wipe down the surface with a dry or lightly damp soft cloth if needed. Store the mat flat or loosely rolled in a dry space, away from direct sunlight or moisture. Hydrate after each session to support lymphatic flow and cellular recovery.

A black and white photograph of a person's legs from the mid-thigh down to the feet. The person is standing on a light-colored surface, with their right leg forward and slightly bent at the knee, and their left leg back and straight. The lighting is dramatic, highlighting the contours of the muscles and the texture of the skin. The background is a plain, light gray.

Muscle Soreness and Recovery

Muscle soreness after an intense workout is caused by microtears in muscle fibers, which lead to inflammation, stiffness, and fatigue. While this process is necessary for muscle adaptation and growth, excessive soreness can slow recovery and impact training consistency.

Support your active lifestyle with effective, non-invasive light therapy.

Muscle soreness is a natural response to intense training. It results from microtears in muscle fibers, causing inflammation, stiffness, and fatigue. While this repair process is essential for growth, prolonged soreness can hinder performance and recovery.

Flowlight LED Mat 1300 Two Waves delivers high-density red and near-infrared light directly to muscles and joints, accelerating recovery and helping you bounce back faster.

Boosting Recovery through Cellular Energy

Red (660nm) and near-infrared (850nm) wavelengths penetrate skin and muscle tissues and are absorbed by the mitochondria—the “power plants” of the cells. This stimulates the production of adenosine triphosphate (ATP), which is critical for:

- Muscle fiber repair
- Tissue regeneration
- Sustained performance during repeated sessions

Why it matters

More ATP means faster recovery between workouts, improved endurance, and a lower risk of overtraining injuries.

Reduced Inflammation and Increased Blood Flow

Near-infrared light (850nm) promotes vasodilation—widening of blood vessels—which increases the flow of oxygen-rich blood and nutrients to sore areas. This helps flush out lactic acid and other metabolic waste more efficiently.

Result

- Faster muscle relaxation
- Less stiffness
- Reduced inflammation post-exercise

Enhanced Microcirculation

Red light improves microcirculation in surface tissue, further supporting nutrient delivery and detoxification. This dual-action on both surface and deeper layers makes the mat particularly effective for large muscle groups like:

- Hamstrings
- Quadriceps
- Glutes
- Lower back


Recommended Use for Muscle Recovery

Positioning

Lay the mat flat and lie directly on it, with the target area in full contact with the surface. The soft design makes it ideal for recovery in a resting, reclined position—whether after a workout or before sleep.

Suggested Use

- Session Duration: 15–30 minutes
For post-training recovery and inflammation reduction
- Short Sessions: 10–15 minutes
For light soreness, warm-up, or daily maintenance
- Frequency: 3–5 times per week
Daily use is recommended for athletes or those in intense training blocks
- Best Time of Day:
 - Post-workout – to speed up recovery
 - Evening or before bed – to relax muscles and prepare for overnight repair

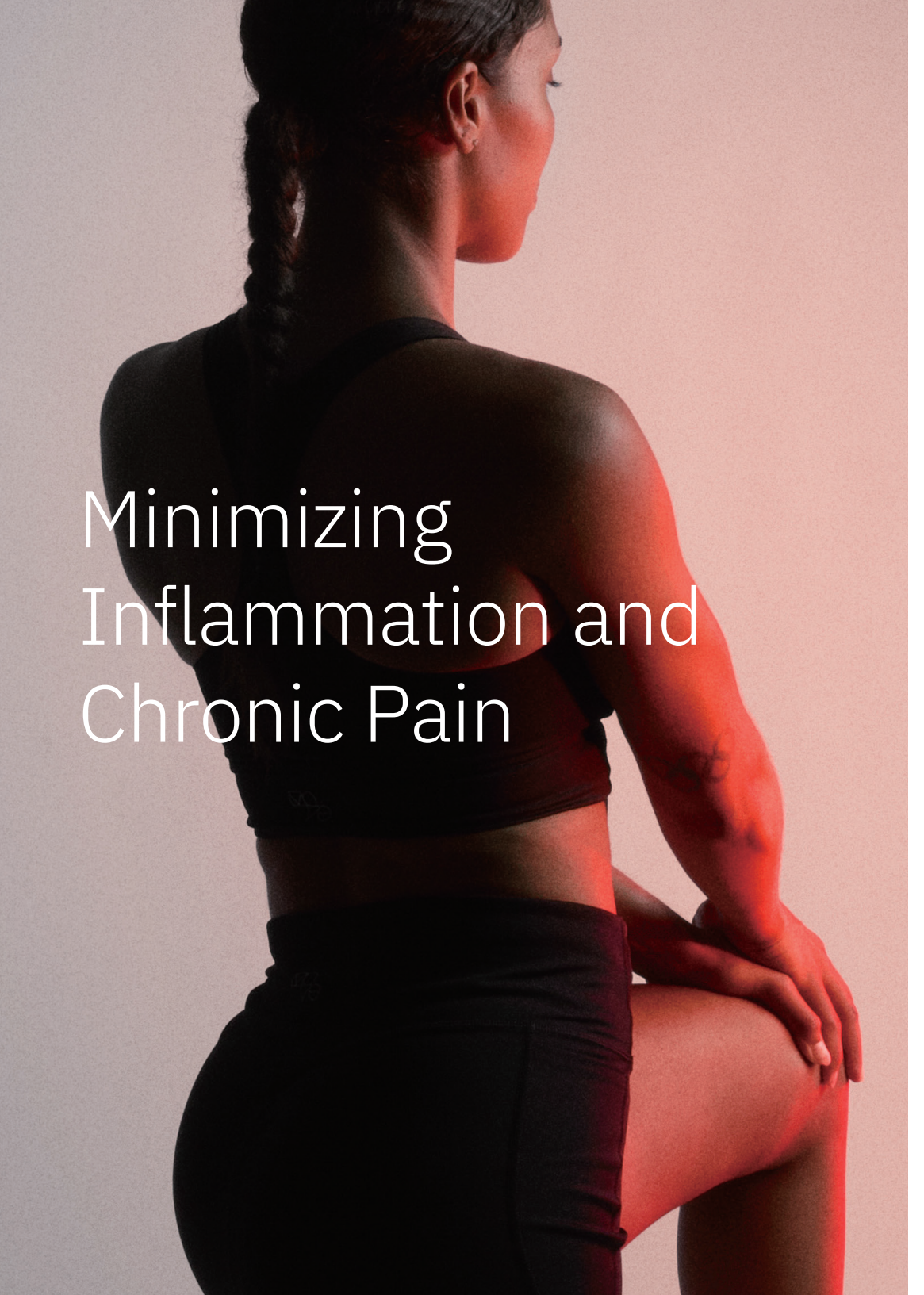Minimizing Inflammation and Chronic Pain

Restore mobility, reduce pressure, and calm the nervous system.

Joint pain. Chronic tightness. Recurring discomfort. These issues often stem from low-grade, unresolved inflammation—deep in connective tissue, fascia, and joints. With regular use, Flowlight LED Mat 1300 Two Waves helps reduce these symptoms by improving circulation, modulating inflammation, and stimulating deep tissue repair.

Unlike passive rest or pain-masking treatments, light therapy targets inflammation at its biological root.

Reducing Inflammation at the Source

Chronic pain is often sustained by persistent inflammation in muscles, joints, and connective tissue. When cells are under stress, they release pro-inflammatory signals that restrict mobility and slow down healing.

Near-infrared light (850nm) penetrates several centimeters into the tissue, where it:

- Increases ATP production
- Reduces oxidative stress
- Modulates inflammatory signaling pathways
- Stimulates repair in joint capsules and tendons

Improving Circulation for Deeper Healing

By expanding capillaries and enhancing blood flow, both red and NIR light improve oxygen delivery to damaged tissue and help remove built-up waste and inflammatory byproducts.

- Red light (660nm): improves circulation in skin and surface-level tissue
- NIR (850nm): restores perfusion in fascia, ligaments, and deep muscle layers

Result

- Reduced joint stiffness
- Decreased tension in overused areas
- Better range of motion and movement quality

Nervous System Regulation

Chronic pain often involves the central nervous system. Regular light therapy has been shown to reduce sympathetic overdrive (“fight or flight”) and support parasympathetic tone (“rest and repair”). This is especially true with full-body formats like the Flowlight Mat, which encourages whole-system relaxation during use.

Why this matters

Pain is not just physical—it’s also neurological. Light therapy can help retrain pain perception, reduce central sensitization, and break long-term tension patterns.

Recommended Use for Inflammation and Pain Relief

Positioning

Lie on the mat with the affected joints or tension zones in direct contact. If needed, support joints with towels or pillows to optimize alignment.

Suggested Use

- Mode: IR/NIR (660 + 850nm)
- Intensity: 100% (adjust if overly sensitive)
- Duration: 20–30 minutes
- Frequency: Daily for chronic conditions, or 3–5x per week for flare-ups
- Best time of day: Evening or post-activity, when inflammation tends to peak

Common Target Areas

Condition	Recommended Setup
Neck & shoulder tension	Upper back flat, arms slightly elevated
Lower back stiffness	Spine centered, knees bent for pelvic tilt
Hip or SI-joint pain	Rotate gently to increase contact on one side
Knee discomfort	Back of knee pressed gently into mat
Elbow/wrist issues	Entire arm placed flat, palm up

Regulating Circadian Rhythm and Combating Seasonal Affective Disorder (SAD) with Light Therapy

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Regulating Circadian Rhythm and Combating SAD with Light Therapy

Light is not just for recovery. It's your body's most powerful clock.

The body's internal rhythm—known as the circadian cycle—is regulated by light exposure. When this rhythm is disrupted, sleep suffers, mood declines, and hormonal balance shifts. Flowlight LED Mat 1300 Two Waves delivers wavelengths that help reset this cycle, especially during darker months or irregular schedules.

Red and near-infrared light help reinforce natural energy patterns, promote melatonin production, and counteract seasonal mood dips.

Synchronizing the Body's Internal Clock

Circadian rhythms are controlled by the suprachiasmatic nucleus (SCN), which responds to light signals from the body. While blue light (such as sunlight or screens) directly stimulates wakefulness, red and near-infrared light help stabilize the circadian rhythm without overstimulating the brain or suppressing melatonin.

Morning sessions using red light signal alertness and support healthy cortisol rhythms.

Evening sessions using near-infrared light promote cellular relaxation and support melatonin release for deeper sleep.

Combating Seasonal Affective Disorder (SAD)

SAD is a form of depression triggered by insufficient light exposure, typically in winter. This results in reduced serotonin levels, disrupted sleep, and lower energy. Unlike artificial blue-light lamps, Flowlight LED Mat 1300 Two Waves delivers biologically calming wavelengths that support mitochondrial health, mood regulation, and sleep quality.

- Red light (660nm) may improve alertness and energy without overstimulation.
- Near-infrared (850nm) influences deeper neuroendocrine systems, helping regulate mood-related neurotransmitters.

Melatonin, Sleep Quality & Hormonal Support

Evening use of near-infrared light has been shown to support natural melatonin production by reducing cellular stress and shifting the nervous system into a parasympathetic state. This leads to:

- Easier sleep onset
- Fewer nighttime awakenings
- Deeper and more restorative sleep cycles

Why this matters

Circadian misalignment is linked not just to poor sleep, but also to impaired recovery, mood swings, and immune dysregulation. Supporting light-driven hormonal patterns is essential for long-term wellness.

Recommended Use for Sleep & Energy Regulation

Morning Sessions (Red light only)


- Mode: RED or RED/NIR
- Time: Within 1 hour of waking
- Duration: 10–15 minutes
- Positioning: Lie on the mat with your back exposed, arms by your side
- Goal: Increase alertness and support cortisol rhythm

Evening Sessions (NIR only):

- Mode: NIR
- Time: 30–90 minutes before bedtime
- Duration: 15–20 minutes
- Positioning: Full-body contact
- Goal: Promote melatonin release and deeper rest

Frequency:

- Daily use is recommended, especially during the winter months or in shift work.

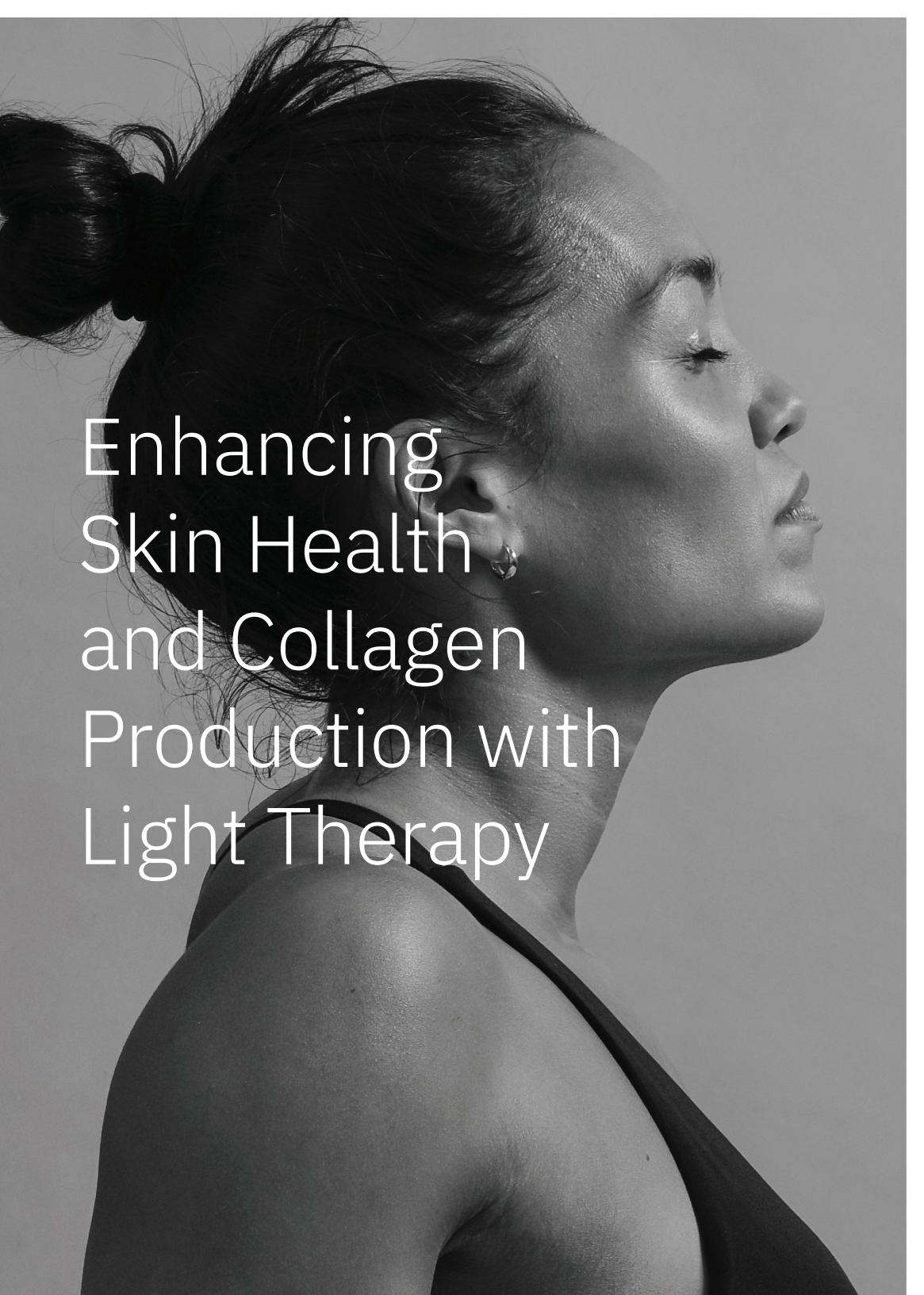A black and white photograph of a woman's profile, facing right. Her hair is dark and pulled back into a bun. She is wearing a dark, possibly black, top. The lighting is soft, highlighting the contours of her face and neck. The background is a plain, light gray. Overlaid on the left side of the image is the text "Enhancing Skin Health and Collagen Production with Light Therapy" in a white, sans-serif font.

Enhancing Skin Health and Collagen Production with Light Therapy

Healthy skin starts below the surface.

Skin aging, dryness, uneven tone, and scars are not just surface issues – they are cellular responses to oxidative stress, reduced collagen production, and impaired microcirculation. Flowlight LED Mat 1300 Two Waves delivers targeted 660nm red light to stimulate repair and regeneration in the skin's upper layers.

Regular use improves texture, elasticity, and tone – naturally and non-invasively.

Boosting Collagen and Elastin Production

Collagen is the structural protein responsible for firm, resilient skin. As we age or experience stress, collagen production slows, leading to wrinkles, sagging, and dullness.

Red light at 660nm is absorbed by fibroblast cells in the dermis, stimulating increased production of:

- Type I and III collagen
- Elastin fibers
- Growth factors that support healthy cell turnover

Result

- Firmer, smoother skin
- Reduced appearance of fine lines
- Improved elasticity over time

Accelerating Healing and Reducing Scars

Whether from acne, injury, or surgical healing, red light therapy supports tissue regeneration by:

- Enhancing local blood flow
- Reducing oxidative stress
- Supporting healthy inflammation resolution

Near-infrared light (850nm) penetrates deeper, accelerating regeneration in underlying tissues and promoting new cell formation beneath scar tissue.

Soothing Inflammation and Redness

Chronic redness, rosacea, or post-inflammatory hyperpigmentation are often driven by ongoing irritation in superficial blood vessels and skin cells.

Red and near-infrared light reduce inflammatory cytokine activity and improve skin oxygenation. This calms the tissue and supports clearer, more even skin tone.

Enhancing Blood Flow for a Radiant Complexion

Both red and near-infrared light stimulate angiogenesis – the formation of new capillaries. This improves nutrient delivery to the skin, removes cellular waste, and results in a more vibrant, healthy glow.

Recommended Use for Skin Rejuvenation

Treatment Guidelines:

- Mode: RED or IR/NIR (for deeper support)
- Intensity: 80–100% depending on sensitivity
- Duration: 10–20 minutes per session
- Frequency: 5–7 times per week for visible changes
- Best time of day: Morning for revitalization, evening for repair

Positioning Tips:


Lie flat on the mat with full-skin contact. For targeted face or décolleté treatment, position the upper body on the mat and adjust arm placement to open the chest and neck area. Ensure skin is clean and free of oils or lotion.

What to Expect Over Time

Week	Visible Changes
1-2	Increased hydration, glow, and skin softness
3-4	Reduced redness, more even tone
4-6	Noticeable firmness and improved elasticity
6+	Fewer fine lines, improved texture, long-term resilience

Note


Light therapy works with the skin’s natural renewal cycle. Consistency and patience are key. Most users see visible improvement within 4–6 weeks of regular use.

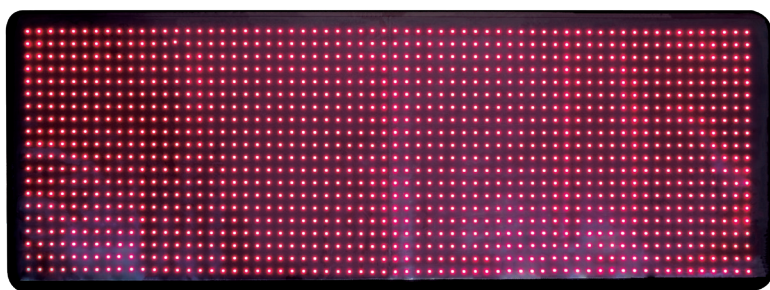hello@flowlife.com
www.flowlife.com
Runebergsgatan 7a
113 51 Stockholm

Manual

Flowlight LED Mat 1300 Two Waves

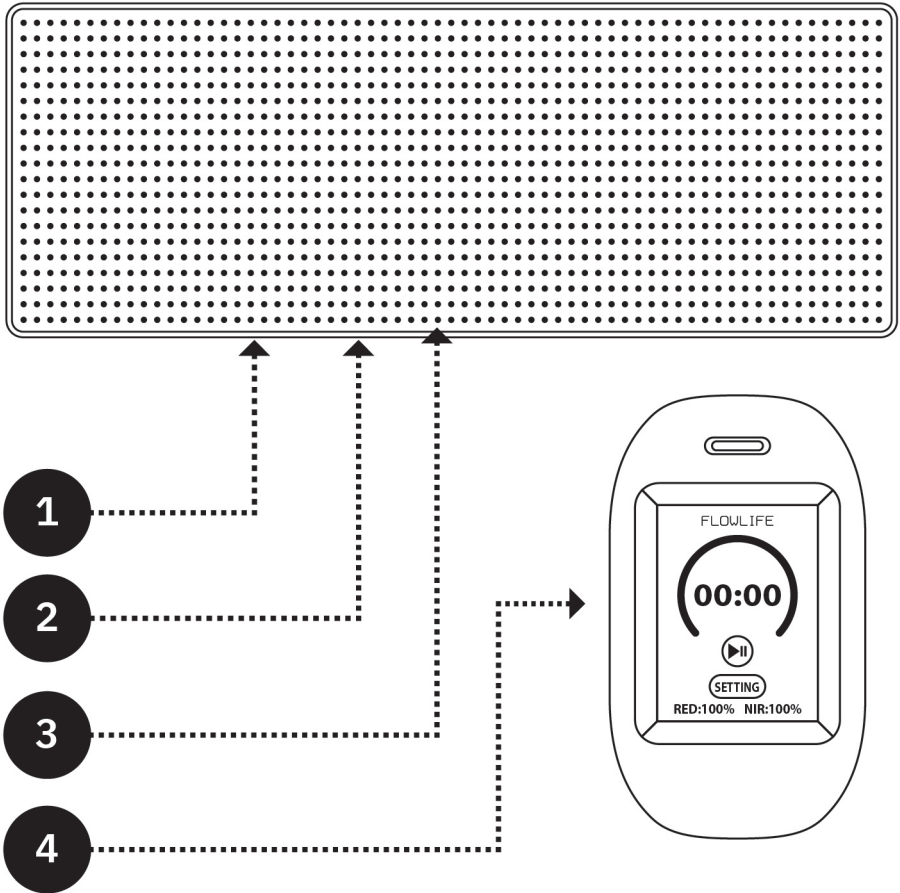1. Fabric material - can be brushed but no water
2. PC cover on the surface
3. Red lights and NIR lights
4. Remote control with LCD screen

ENGLISH

INTRODUCTION

Welcome to your Flowlight Mat user manual. This guide will help you get started with your new red and near-infrared light therapy device and ensure you understand how to use it safely and effectively.

In this document, Flowlight Mat is referred to as the product, mat or Flowlight.

The following safety instructions must be used for the product. This reduces the risk of fire, electric shock, property damage, and personal injury. Read all safety instructions before using the Flowlight.

SAFETY INSTRUCTIONS

SAFETY HAZARD – REDUCE RISK OF ELECTRIC SHOCK

- Never have the product plugged into an electrical outlet or other power source when cleaning the product.
- Never use the product near water, such as in a bath, shower, or similar.
- Never place the product anywhere where it risks being pulled into the water.
- Always store the product in a dry place, without contact with water or other liquids.
- Immediately unplug the product if it comes in contact with water.

DANGER – REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK AND PERSONAL INJURY

- Never leave the product unattended after connecting it to a power source.
- The product must never be covered with fabric or similar materials during use. This can lead to overheating, fire, electric shock, and personal injury.
- Provided the user has been properly informed about safe usage and the risks outlined in the safety instructions, children aged 15 and above, people with various disabilities, and those with limited knowledge or experience may use the product.
- Never use the product for any purpose other than what is described in this safety instruction. It should only be used for light therapy on close distance.
- Never use product accessories from any manufacturer other than the original.
- Never use the product if there is apparent or suspected damage to the device, charging cord, plug, or power source.
- If the product has been in contact with water, has been dropped, or damaged in any other way, always contact customer service for inspection and potential repair.
- Only use the product once you have read the customer service information and recommendations.
- Never use the power cord as a handle to carry the product.
- Never store the product or cable on hot surfaces.
- Never insert anything into the product openings.
- Never use the product where aerosol products are used or where oxygen is administered.
- Children may not use, clean, or maintain the product without adult supervision.
- Never attempt to repair a broken product yourself. Never open or unscrew it. Customer service must always be contacted in case of apparent or suspected damage.
- Always consult your doctor before using the product if you need any clarification.

RISK GROUPS

If you belong to any of the following risk groups, ALWAYS consult your doctor before using Flowlight. Use the product only with your doctor's approval.

- Recent Burns: Avoid using red light therapy on recent burns.
- Heat Sensitivity: Red light therapy generates heat, which may pose risks in some conditions.
- Light Sensitivity: Individuals with extreme photosensitivity or who have previously reacted to red light therapy should avoid its use.
- Malignant Cancers: There is a theoretical risk that red light therapy may increase cell proliferation in malignant cancers.
- Hyperthyroidism: Avoid applying red light therapy on the neck and upper chest.
- Pregnancy: Avoid using red light therapy on the abdomen during pregnancy. Consult with a healthcare provider first.
- Eye Disease: Consult with a healthcare provider if you have any eye condition before using red light therapy.
- Fever or Infection: Red light therapy may worsen hyperthermia caused by fever or infection.
- Systemic Lupus Erythematosus (SLE): Increased risk of adverse reactions.
- Severe Bleeding or Blood Loss: Some studies suggest red light therapy may slow blood clotting.
- Photosensitizing Medications: Consult with a healthcare provider if taking medications that increase sensitivity to light.
- Melasma / Hyperpigmentation: Start slowly if concerned about pigmentation changes.
- Safety Glasses: These are provided and must be worn when facing the lights.
- Detox Symptoms: May cause detox-related symptoms initially. Start with shorter sessions and increase gradually.
- Pre-existing Conditions / Medications: Consult with your doctor if you have any health concerns or use photosensitizing medication.

RECYCLING AND SAFETY

- Connect and remove the cord carefully.
- Flowlight must be recycled at an approved electronic waste facility.
- Avoid exposing Flowlight to extreme conditions or temperatures above 50°C to prevent fire hazards.
- Never disassemble Flowlight.

WARRANTY

If your product breaks within two years from the date of purchase due to material or manufacturing defects, Flowlife will repair or replace necessary components or the product free of charge. Flowlife AB is not responsible for warranty or return cases for products purchased from third-party retailers.

The warranty is not valid if:

- the product is damaged due to improper handling, misuse, accident, or transport
- the product has been subject to unauthorized repair
- the product has not been used according to Flowlife's instructions
- the damage exceeds the product's value
- the product deteriorated due to improper storage
- dated proof of purchase cannot be provided

INSTRUCTIONS FOR USE AND MAINTENANCE

Clean the light by unplugging it and wiping the surface with a cloth or paper towel. Avoid chemicals or bleach. Ensure the light is completely dry before use. Indoor use only.

NIR LEDs: Near-infrared light is invisible to the naked eye. Use a camera to check whether the LEDs are functioning.

INSTRUCTIONS FOR SETUP, USE AND STORAGE

To use the Flowlight, follow the instructions below.

FLOWLIGHT PANELS SERIES SETUP

In the package: Flowlight LED Mat 1300 Two Waves LED Mat. 1x UL-certified power driver (24V5A). LCD display controller. 1x Protective glasses. Instruction manual. Travel carry bag.

SET IT UP

Use the mat on a firm, flat surface such as a massage table, yoga mat, or floor. Your skin should be clean, dry, and free from oils or lotions, as these can reduce light penetration. Remove thick clothing from the treatment area to allow for direct light exposure. If you're treating near the face, use the included protective glasses.

STARTING YOUR FLOWLIGHT MAT

Connect the male connector of the remote cable to the cable attached to the Flowlight Mat. Make sure the slots on both connectors are aligned, then push them together firmly and secure the connection by tightening the locking ring. Connect the female connector of the remote cable to the power adapter cable. Align the slots, press firmly until fully connected, and tighten the locking ring.

Connect Flowlight to a power source and switch it on using the LCD remote.

PROGRAMMING

Control all functions through the LCD screen on the remote. The pause button allows you to start or stop the treatment. When the session begins, the timer counts down.

ENDING YOUR SESSION

Turn off the device by turning it off with the remote. Unplug the device. Allow the mat to cool. Wipe clean with a soft, dry cloth and store safely.

REMOTE CONTROL

The remote control can also control all settings.

Setting the timer: 1–30 minutes

Adjusting light intensity: 1%–100%

Pulse setting: 1–10,000 Hz

STANDARD USE

First-Time Users: Begin with 10-minute sessions per area.

General Use: 10–30 minutes at 15–30 cm distance.

Allow 6 hours between treatments on the same area.

Do not exceed 30 minutes per session.

Use 4–12 times weekly with at least one rest day.

Wear protective glasses when facing the light.

DURING YOUR SESSION

Sit back, relax, and let the Flowlight work.

Adjust settings as needed for comfort and effectiveness.

STORAGE

Store in a cool, dry place.

Unplug when not in use.

2 WAVELENGTHS

RED/NIR: 660 nm + 850 nm

RED: 660 nm only NIR: 850 nm only

TECHNICAL SPECIFICATIONS

Spectrum: 660nm (Red) & 850nm (NIR)

Light Intensity: 1%–100%

Input Voltage: 100V–240V,

Frequency: 50–60Hz

Rated Power: 1280 W

Actual Power: approx. 97 W

Dimensions: 1575X600mm

Timer: 1 - 30 minutes

Gross Weight: 2,95 kg per unit

Color: Black

Shell Material: Fabric

Irradiance: approx. 130 mW/cm² at surface level

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